

Optical Imaging Data

Labeo Technologies, September 2020

Recommended Hardware

- Any recent CPU should work without any issues.
- 32 GB of RAM.
- SSD Hard drive.
- The analysis software uses Matlab and toolboxes: wavelength & signal processing.

How raw data is organized:

- Each recordings have their own folder [A]
- In this folder, there will be 3 different type of files:
 - **Info.txt**: contains all the parameters of the system for that specific recording
 - **ai_#####.bin**: binary file containing analog inputs
 - **img_#####.bin**: binary file containing images

Name	Date modified	Type	Size
info	5/21/2020 3:54 PM	Text Document	1 KB
ai_00000.bin	1/31/2019 11:48 AM	BIN File	157,266 KB
img_00021.bin	1/31/2019 11:48 AM	BIN File	202,755 KB
img_00020.bin	1/31/2019 11:47 AM	BIN File	524,295 KB
img_00019.bin	1/31/2019 11:47 AM	BIN File	524,295 KB
img_00018.bin	1/31/2019 11:47 AM	BIN File	524,295 KB
img_00017.bin	1/31/2019 11:47 AM	BIN File	524,295 KB
img_00016.bin	1/31/2019 11:47 AM	BIN File	524,295 KB
img_00015.bin	1/31/2019 11:47 AM	BIN File	524,295 KB
img_00014.bin	1/31/2019 11:47 AM	BIN File	524,295 KB
img_00013.bin	1/31/2019 11:46 AM	BIN File	524,295 KB
img_00012.bin	1/31/2019 11:46 AM	BIN File	524,295 KB
img_00011.bin	1/31/2019 11:46 AM	BIN File	524,295 KB
img_00010.bin	1/31/2019 11:46 AM	BIN File	524,295 KB
img_00009.bin	1/31/2019 11:46 AM	BIN File	524,295 KB
img_00008.bin	1/31/2019 11:46 AM	BIN File	524,295 KB
img_00007.bin	1/31/2019 11:46 AM	BIN File	524,295 KB
img_00006.bin	1/31/2019 11:45 AM	BIN File	524,295 KB
img_00005.bin	1/31/2019 11:45 AM	BIN File	524,295 KB
img_00004.bin	1/31/2019 11:45 AM	BIN File	524,295 KB
img_00003.bin	1/31/2019 11:45 AM	BIN File	524,295 KB
img_00002.bin	1/31/2019 11:45 AM	BIN File	524,295 KB
img_00001.bin	1/31/2019 11:45 AM	BIN File	524,295 KB
img_00000.bin	1/31/2019 11:45 AM	BIN File	524,295 KB

Update Notes:


- The names of the files generated by ImagesClassification changed in a recent update. Be aware of this when using the following scripts, there might be errors due to this change.
 - Old file names: rChan.dat, gChan.dat, yChan.dat, fChan.dat and sChan.dat
 - New file names: red.dat, green.dat, yellow.dat, fluo_xx.dat and speckle.dat

Opening data from Matlab


- An open source library with the basic opening/filtering functions is supplied by Labeotech.
- This library can:
 - Open and classify the different channels
 - Compute HbO and HbR from the non-fluorescent channels (*to be released soon*)
 - Do the hemodynamic correction on fluorescence data
 - Compute $\Delta F/F$ for fluorescence channels
 - Compute blood flow from speckle acquisitions

Data organization in .bin files:

- img files:
 - Channels are interlaced during recordings.
 - Each file contains X images (depending on the size of the region of interest, this number can change. Files will always be around 500 MB)
 - Once the first file is full, a new img_ file is created, until the end of the acquisition.

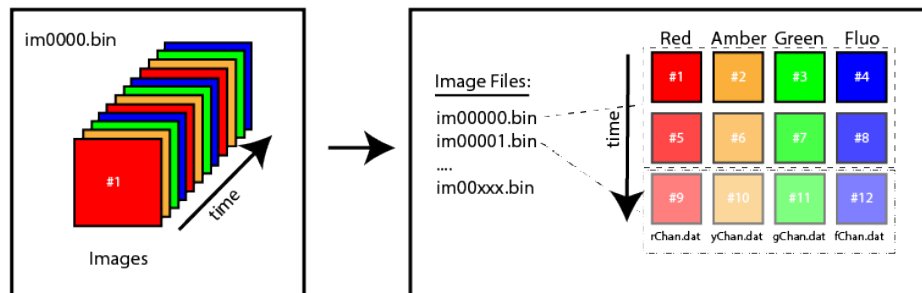

- ai files:
 - The data pattern on those file is a bit more complex.
 - Analog inputs are saved once every seconds (and at the end of the acquisition)
 - The transmission of the data from each AI channel is interlaced when received from the acquisition board.
 - For each transmission, the new chunk of data is stacked at the end of the previous one. Once the file size limit is reach, a new ai_####.bin file is created.

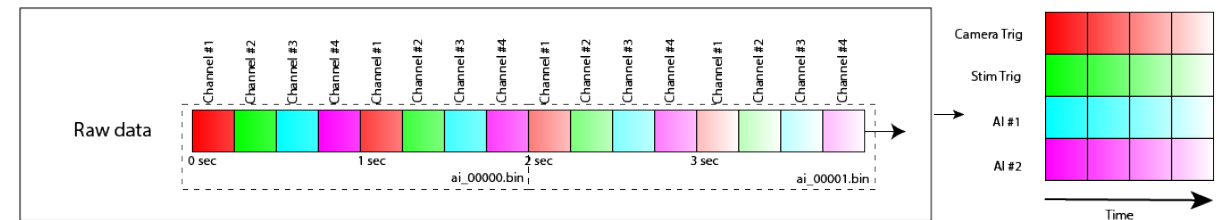

Pre-Analysis (to separate channels):





Fonction to put each recorded channel in its own file:

`ImagesClassification(sDataFolder, sSaveFolder, iSpatialBinning, iTemporalBinning, blgnoreStim, bROI);`

1. First parameter is the path to the folder containing data to be analysed;
2. Second is the path to where to save the results.
3. Third parameter is for the spatial binning option (set to 2 for a 2x2 binning; 1 for no binning or 1x1.);
4. Fourth parameter is for temporal binning (through averaging). Set to 1 for no binning, and to any integer number for binning (i.e.: 3 for a 3 frames averaging)
5. Fifth parameter is to tell the software to consider the acquisition as a Resting state protocol (meaning without stimulation). (set to 1 to ignore stimulation signal; set to 0 to do as required based on the stimulation signal).
6. Sixth parameter is to be able to set a post-acquisition Region of Interest. It can help to reduce the size of an acquisition. Set to 0 to ignore.

Normalization filtering

- Let's suppose a temporal signal from one pixel (figA)
- The filter uses two frequencies to normalize the signal:
 - A low frequency filter
 - A high frequency filter
- In the case of the pixel presented before, we would get these two signals (figB)
- The normalized output would then be as presented in (figC)

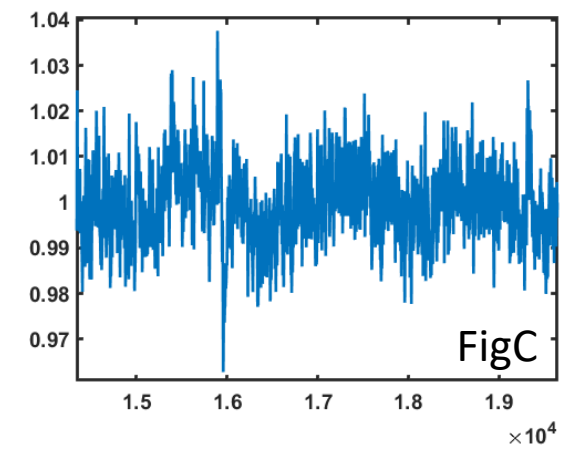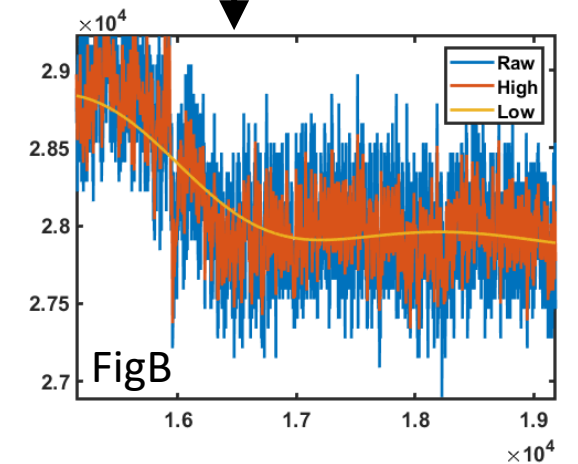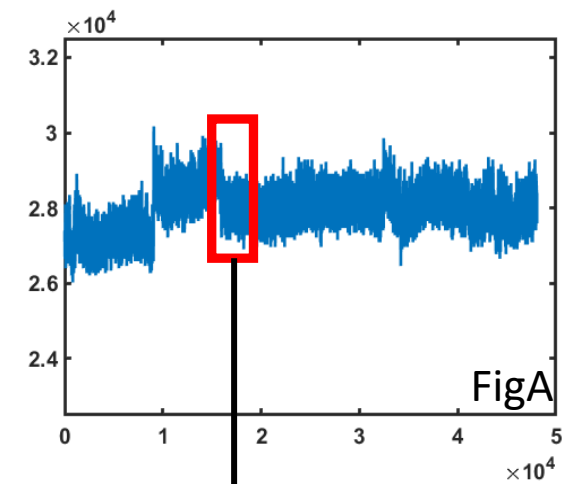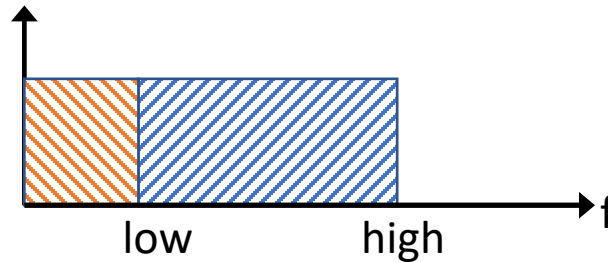

Normalization:

(Must be runned after ImagesClassification.m)

To do the intrinsic channel normalisation (see following slide for details)

NormalisationFiltering(*FolderData, FileData, lowFreq, highFreq, bDivide*)

1. Path of the folder containing the data file to be normalized.
2. Can be either: the name of the file to be normalized or the matrix containing the data.
3. Low frequency cut-off
4. High frequency cut-off
5. If set to 1: the returned data will be of the form: $Data_{HF} / Data_{LF}$

Otherwise: $Data_{HF} - Data_{LF}$

Ex: `dat = NormalisationFiltering(pwd, 'red', 1/120, 1, 1);`

Ex2: `dat = NormalisationFiltering(pwd, dat, 0, 3, 0);`

Hemodynamic correction

(Must be runned after ImagesClassification.m)

To do the hemodynamic correction on Fluorescence

HemoCorrection(*Folder, varargin*):

- 1- Folder: Folder containing the dataset to work with.
- 2- Varargin a)-> if empty: a dialog box will be prompt to ask user which channels to use to do the correction.
b) -> cell array of string: to specify which channels to use. (Red, Green and/or Yellow)

Ex: HemoCorrection(pwd, {'Red', 'Green'});

To do the fluorescence normalisation:

Dat = NormalisationFilter(*FolderData, FileData, lowFreq, highFreq, bDivide*)

- See slides 8 and 9 for explanations.
- Ex: dat = NormalisationFilter(pwd, 'fluo_475', 0.3, 3, 1);

How to run analysis (for Blood Flow):

(Must be runned after ImagesClassification.m)

To get a quantitative measure of blood flow

Ana_Speckle(*sFolderName*, [])

1. First parameter is the path to the folder containing data;
2. Second parameter is not used at this time, must be set to [];

To get a mapping of stdev of speckle images

SpeckleMapping(*folderpath*, *sType*)

1. First parameter is the path to the folder containing data;
2. Second parameter is to set which type of stdev to be computed: « Spatial » or « Temporal »

How to see results:

- Manually:

- In Matlab (change the working directory to where the data is stored), then:

```
Info = matfile('fluo_475.mat');
```

```
fid = fopen('fluo_475.dat');
```

```
dat = fread(fid, inf, 'single');
```

```
fclose(fid);
```

```
dat = reshape(dat, Info.datSize(1,1), Info.datSize(1,2),[]);
```

```
Imagesc(dat(:,:,1)); %or any other frame #
```

- This can be done for any .dat file:
 - sChan.dat or speckle.dat for speckle
 - rChan.dat or red.dat, gChan.dat or green.dat and yChan.dat or yellow.dat for any of the three colors channels
 - fChan.dat (or fluo_xx.dat) for fluorescence

Example:

- Step 1: Pre-analyse

```
ImagesClassification(pwd, pwd, 1, 1, 0, 0);
```

- Step 2: Hemodynamic correction

```
dat = HemoCorrection(pwd,{'red','green','yellow'});
```

- Step 3: Normalization

```
Dat = NormalisationFiltering(pwd, dat, 0.3, 3,1);
```

- Step 4: Visualisation

```
figure;
```

```
for ind = 1:size(dat,3)
```

```
    imagesc(dat(:, :, ind), [0.95 1.05]);
```

```
    title(int2str(ind));
```

```
    pause(0.1);
```

```
end
```